


BAG - Test Systems for Monitoring Cleaning Efficacy

(available on request only)

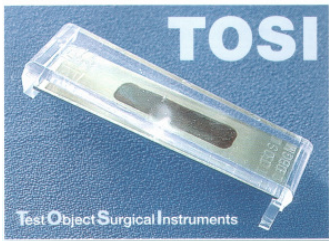


TOSI®

Test device for monitoring cleaning efficacy of washer-disinfectors when reprocessing surgical instruments. Test soil correlates to human blood acc. to ISO 15883.

Order-No.: 8302

Packaging unit: 12 tests/box

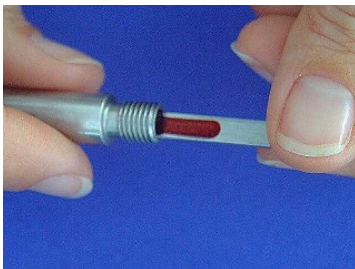

TOSI®-LumCheck

Test device for monitoring cleaning efficacy of washer-disinfectors when reprocessing MIC instruments. Test device simulates hollow loads. Test soil correlates to human blood acc. ISO 15883.

TOSI®-LumCheck

Order-No.: 8303


Packaging unit: 25 tests/box



LumCheck test device

Order-No.: 8305

Packaging unit: 1 device/box

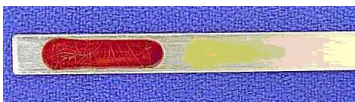

TOSI®-FlexiCheck

Test device for monitoring cleaning efficacy of endoscope washers. Test device simulates endoscope channels. TOSI®-FlexiCheck combines two different test soils on one test strip to check the removal of blood and polysaccharides.

TOSI®-FlexiCheck

Order-No.: 8298


Packaging unit: 25 tests/box



FlexiCheck test device

Order-No.: 8299

Packaging unit: 1 device/box

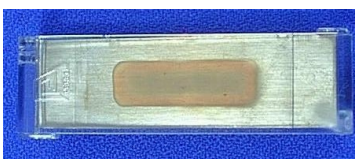

TOSI® - Gold

Test device for monitoring cleaning efficacy of washer-disinfectors when reprocessing surgical instruments with strong protein dissolving chemicals. Test soil correlates to denaturated water insoluble proteins.

Order-No.: 8309

Packaging unit: 12 tests/box



SonoCheck

Test for monitoring the energy in ultrasonic cleaners used for reprocessing of surgical instruments.

Order-No.: 8300

Packaging unit: 30 tests/box



HemoCheck-S

Test kit for qualitative detection of blood residues on surgical instruments and surfaces. Blue colour change within 1 minute with detection limit of 0.1 µg blood. No misinterpretation caused by fingerprints. Longer and smaller swabs for lumen available.

Order-No.: 8295

Packaging unit: 12 tests/box



Pyromol Test

Test kit for qualitative detection of protein residues on surgical instruments and surfaces. Blue colour change within 1 minute with detection limit of 1 µg protein.

Order-No.: 8291

Packaging unit: 20 tests/box



HemoCheck-E / Pyromol-E

Test kit for qualitative detection of blood residues or proteins in biopsy channels. Swab length 2.55 meters. Colour change within 1 minute with detection limit of 0.1 µg blood or 1 µg protein.



HemoCheck-E

Order-No.: 82955

Packaging unit: 6 tests/box

Pyromol-E

Order-No.: 82915

Packaging unit: 6 tests/box



Efferzyme

Enzymatic cleaning tabs for ultrasonic baths. Non-corrosive, will not harm instruments

Order-No.: 8301

Packaging unit: 52 tablets/box